


# oYo-Link®

## Click Chemistry & Chemical Conjugations

Use various chemistries with your antibody of choice. Perform click-reactions prior to antibody conjugation for optimal reaction conditions & higher efficiency.

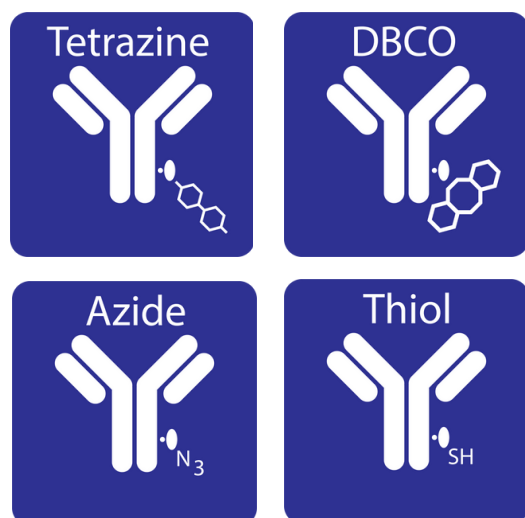

- ✓ Over **90% more efficient** click reactions.
- ✓ **Site-specific, covalent labeling** guarantees no interference with antibody performance.
- ✓ **<30 seconds “hands-on” time, 2 hours total.** Just mix & illuminate.
- ✓ **Compatible with all common buffers,** including BSA, Tris & Azide.

oYo-Link® enables site-specific labeling to the heavy chain of compatible antibodies with up to two click chemistry labels per antibody.

The click reaction occurs prior to adding the antibody, therefore the reaction is not limited by the antibody concentration, resulting in a highly efficient reaction.

Label as little as 1 µg of antibody at a time. Work with concentrations as low as 50 µg/mL.

